# BECKEN

[you can trust]

## ***Baby and Toddler Scale***

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MANUAL DE INSTRUCCIONES

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Dear client,

Thank you for purchasing this product.

The quality of this appliance has been certified by several laboratory tests. This instructions manual is provided to simplify the use of the appliance as well as to guarantee a safe use.

Please read carefully the instructions and safety measures described in this manual.

We sincerely hope you will be satisfied when using your new appliance.

## INDEX

<b>1. SAFETY INSTRUCTIONS</b>	<b>4</b>
1.1. Humidity and Water	5
1.2. Maintenance & Cleaning	5
1.3. General Precautions	6
<b>2. OPERATING INSTRUCTIONS</b>	<b>7</b>
2.1. Product Description	7
2.2. Using the Product	7
<b>3. TECHNICAL SPECIFICATIONS</b>	<b>10</b>
<b>4. POST-SALE SERVICE</b>	<b>10</b>
<b>5. ENVIRONMENTAL PROTECTION</b>	<b>11</b>

# BECKEN

## 1. SAFETY INSTRUCTIONS

	<b>WARNING</b>	
	RISK OF ELECTRIC SHOCK DO NOT OPEN	
<b>Warning:</b> to avoid the risk of electric shock, do not remove the cover. The interior does not contain any parts that can be repaired by the user. Ask for assistance from qualified personnel.		

### **Important Warnings**

Before using your appliance read this instructions manual carefully and keep it in a safe place as you may need to consult it in the future.

Before turning on your appliance, and to ensure a safe and correct operation, do not use if:

- It has fallen to the ground;
- You notice any evidence of damage;
- Any anomaly occurs during its operation.

To avoid hazards or deterioration of your appliance through incorrect use it is advisable to read these instructions thoroughly.

Your appliance may not be used for other purposes than the ones intended and solely for household use. Any damage resulting from using the appliance outside this scope, any incorrect use, as well as any manual modifications to the product will not be covered and automatically voids the warranty.

**Assistance**

In case of malfunction, take your appliance to a qualified technician for assistance.

This appliance can be used by children from the age of 8 and by people with reduced physical, sensorial or mental capabilities and/or lack of experience and knowledge if they have been given instructions on using the appliance safely and if they understand the hazards involved. Children may not play with the appliance. Cleaning and maintenance must not be carried out by children without supervision.

These instructions are also available at [www.suporteworten.pt](http://www.suporteworten.pt)

**1.1. Humidity and Water**

Do not use your appliance in humid places. Do not allow the appliance to become wet under any circumstances as this may be dangerous. To avoid fire hazard or electric shock do not expose your appliance to humidity or water. Furthermore, do not place your appliance under water (e.g. for cleaning).

**1.2. Maintenance & Cleaning**

Before cleaning, be sure to switch off the appliance and remove the plug from the main power supply. Always allow the unit to cool down before cleaning.

To clean the exterior of the product use a dry soft cloth only.

- Clean the scale only with water, a mild soap and a soft cloth. Prevent water from getting inside the scale.
- Do not immerse it in water. Avoid abrasive cleaners and solvents.
- When the scale is not used for a long time, please make sure you remove the batteries.

# BECKEN


## 1.3. General Precautions

- Please follow the instructions in the user manual which help to properly install and operate this product and enjoy its advanced features. Please keep the instructions for future reference.
- Read all instructions before using the product.
- Maintenance and service operations are only to be carried out by authorized professionals.
- Any damage caused by manual modifications to the product is not covered by the warranty.
- This scale is designed for weighting Babies and Toddlers.
- Do not use it for commercial purpose.
- Place the scale on a hard and flat surface.
- Do not press or put any object on the platform when not in use.
- Extreme electromagnetic conditions may adversely affect the performance of the scale.
- Keep the scale at a safe distance from such equipment or turn off the equipment when using the scale.
- Once the scale has interference, shut it off and turn it on again; the scale can once again be used for its intended purpose.

## 2. OPERATING INSTRUCTIONS

### 2.1. Product Description

#### Function Key
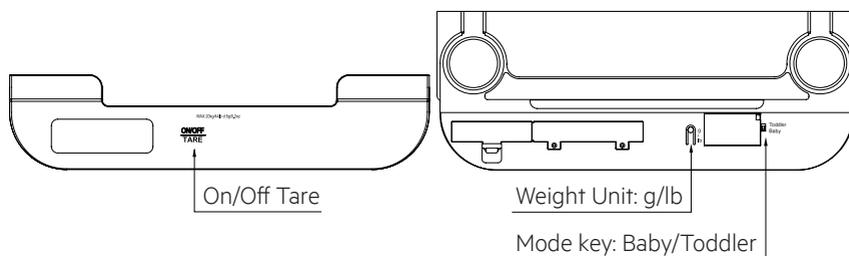
### 2.2. Using the Product

#### Weighing

- Insert batteries, three (3) AA batteries required.
- Press ON/OFF button to turn the power on, all digits will show on the LCD for 5 secs., then "0" g will be shown. The scale is ready for weighing. Keep pressing ON/OFF for 3 seconds to turn off the scale.


- Before weighting, select the mode switch - Baby/Toddler (mode switch is on the opposite side of the LCD panel).
- Before weighting, select weighting unit - g/lb (weighting unit button is on the opposite of the LCD panel).

# BECKEN

## Baby Mode

- If you choose the Baby mode, the word “Baby” will show on the LCD panel, at the top right corner.
- If you choose the Baby mode, please place the scale tray on the scale.

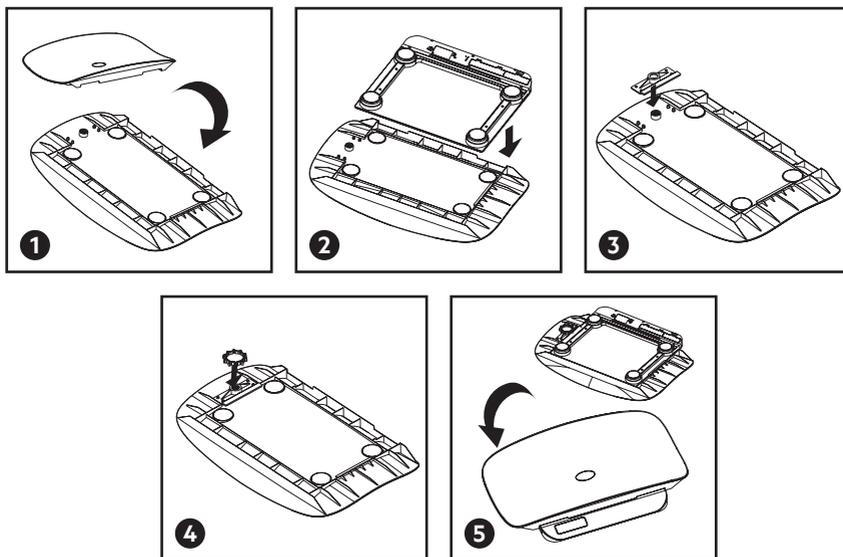- Press ON/OFF/TARE button to turn the power on, all digits will show on the LCD for 5 seconds, then “0”g will be shown. The scale is ready for weighing.

- Carefully place the Baby/Toddler on the scale. The weight value will flash 3 times and then stop.
- Keep pressing ON/OFF/TARE button for 3 seconds to turn off the scale. Otherwise, the power will be shut off automatically within 30 seconds if the LCD=0 and within 20 seconds if the LCD≠0.

### Toddler Mode

- If you choose the “Toddler” mode, press ON/OFF/TARE button to turn on the power. All digits will show on the LCD for 5 seconds, then “0”g will be shown. The scale is ready for weighing.
- Carefully place the Baby/Toddler on the scale. The weight value will show. Remove the Baby/Toddler from the scale. The value will automatically calibrate to “0”g and the scale will be ready for another weighing.
- Keep pressing ON/OFF/TARE button for 3 seconds to turn off the scale. Otherwise, the power will be shut off automatically within 4 minutes if the LCD=0 and within 8 minutes if the LCD≠0.

### TARE Weight Function

“TARE Weight” function allows you to weight several objects separately in the same container. Please check operation procedure below:

- Weight the first object. Its weight will be displayed as mentioned above.
- Press ‘TARE’ button, the display will be reset to zero.
- Add the next object. Please note that only the weight of the second object will be displayed.
- Repeat the procedure for all remaining objects.

### Notice:

- **Under “TARE Weight” function, the total weight cannot exceed the capacity of the scale.**
- **If objects are taken out after setting the “TARE” mode, the LCD will show the weight of the removed objects (in minus). Press ‘TARE’ button to reset to zero.**

# BECKEN

## Zero Calibration

When the scale is turned on the display cannot return to 0g when the object is picked up from the tray. Press ON/OFF TARE button to reset to zero.

## 3. TECHNICAL SPECIFICATIONS

### Function Mode

### Baby Scale/Toddler Scale

Capacity:	Max: 20 Kg / 44 lb < 10 kgs, "G" display >= 10 kgs, "KG" display Graduation: 5 g / 0.2 oz
Power source:	4.5 V (AA x 3)
LCD Panel/view area:	Back light LCD / 26 x 79 (mm)
Overloaded Indication:	-
Low Power Indication:	LO
Tray Dimensions:	500 x 304.5 x 57.6 (mm)
Scale Dimensions:	303 x 340 x 28 (mm)
Glass Plate Dimensions:	250 x 345 x 4 (mm)

## 4. POST-SALE SERVICE

BECKEN has designed this appliance in order to guarantee maximum reliability. However, if a problem occurs do not attempt to open the appliance - there is a risk of electrical shock. If you are not able to solve the problem, please contact the Customer Support Service of any Worten, Modelo or Continente store in order to be assisted. Worten Equipamentos para o Lar provides a 2-year warranty for manufacturing defects from the date of purchase of the appliance.

**WARNING:** Any attempt to repair the appliance without contacting the manufacturer or seller will void the warranty.

All users must be familiar with these instructions. Knowing what could go wrong can help prevent problems from occurring.

## 5. ENVIRONMENTAL PROTECTION

In order to protect the environment we strive to reduce the volume of packaging used, limiting it to three easily sorted materials: cardboard, paper and plastic. The appliance is manufactured using recyclable materials after disassembled by a specialized company. Please comply with local regulations concerning the recycling of all materials.

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*Worten Equipamentos p/ o Lar S.A.