BECKEN

[you can trust]

Body Analysis Scale

*Báscula de Análisis Corporal
Balança de Análise Corporal*

INSTRUCTIONS MANUAL

MANUAL DE INSTRUCCIONES

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BECKEN

*Importado, produzido e distribuído por:
Importado, producido y distribuido por:
Imported, produced and distributed by:*

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Dear client,

Thank you for purchasing this product. The quality of this appliance has been certified by several laboratory tests. This instructions manual is provided to simplify the use of the appliance as well as to guarantee a safe use.

Please read carefully the instructions and safety measures described in this manual.

We sincerely hope you will be satisfied when using your new appliance.

INDEX

1. SAFETY INSTRUCTIONS	4
1.1. Humidity and Water	5
1.2. Cleaning	5
1.3. General Precautions	5
2. OPERATING INSTRUCTIONS	7
2.1. Product Description	7
2.2. Using the Product	8
2.3. Troubleshooting	13
3. TECHNICAL SPECIFICATIONS	15
4. POST-SALE SERVICE	16
5. ENVIRONMENTAL PROTECTION	16

1. SAFETY INSTRUCTIONS

	WARNING	
	RISK OF ELECTRIC SHOCK DO NOT OPEN	
Warning: to avoid the risk of electric shock, do not remove the cover. The interior does not contain any parts that can be repaired by the user. Ask for assistance from qualified personnel.		

Important Warnings

Before using your appliance read this instructions manual carefully and keep it in a safe place as you may need to consult it in the future.

Before turning on your appliance, and to ensure a safe and correct operation, do not use if:

- It has fallen to the ground;
- You notice any evidence of damage;
- Any anomaly occurs during its operation.

To avoid hazards or deterioration of your appliance through incorrect use it is advisable to read these instructions thoroughly.

Your appliance may not be used for other purposes than the ones intended and solely for household use. Any damage resulting from using the appliance outside this scope, any incorrect use, as well as any manual modifications to the product will not be covered and automatically voids the warranty.

Assistance

In case of malfunction, take your appliance to a qualified technician for assistance.

This appliance can be used by children from the age of 8 and by people with reduced physical, sensorial or mental capabilities and/or lack of experience and knowledge if they have been given instructions on using the appliance safely and if they understand the hazards involved. Children may not play with the appliance. Cleaning and maintenance must not be carried out by children without supervision.

These instructions are also available at www.suporteworten.pt

1.1. Humidity and Water

Do not use your appliance in humid places. Do not allow the appliance to become wet under any circumstances as this may be dangerous. To avoid fire hazard or electric shock do not expose your appliance to humidity or water. Furthermore, do not place your appliance under water (e.g. for cleaning).

1.2. Cleaning

Always allow the unit to cool down before cleaning. To clean the exterior of the product, use a dry soft cloth only.

1.3. General Precautions

- Please follow the instructions in the user manual which help to properly install and operate this product and enjoy its advanced features. Please keep the instructions for future reference.
- Read all instructions before using the product.
- Maintenance and service operations are only to be carried out by authorized professionals.
- Any damage caused by manual modifications to the product is not covered by the warranty.
- If this product is used for other purposes than the ones originally intended, or if it is not correctly connected, the product may become damaged and the warranty will be voided.


BECKEN

- The scale must not be used by anyone with a pacemaker or with any other medical implants. The current passed through the body to measure the impedance of its different tissues will affect the proper functioning of such implants.
- The scale is for domestic use only, are not a medical device and should not be used as such.
- Do not use the scale on a wet slippery surface. Do not stand on at the edge of scale to avoid falling over.
- Disabled persons, pregnant women and the elderly are advised to use the scale with the assistance of another person.
- Body fat measurements will vary with the amount of body water and will be affected by the following factors: diet, alcohol intake, exercise, certain prescription drugs or diuretics, illness or, in the case of women, their menstrual cycle.
- Please take off your socks and keep your bare feet clean and dry.
- Advice to use the scale without any clothes or jewellery. This will influence the accuracy of the reading.
- Overeating and water loss will have a negative effect on the measurement.
- Heavy exercise will a have negative effect on the reading.
- Unstable or incorrect measurements can be caused by poor contact between your feet and the metal plate on the scale, please clean any dirt on the scale and follow the provided instructions.
- Keep your balance while weighing, any significant level of movement will cause the measurement to be inaccurate.
- The scale is not guaranteed to work accurately for the following groups:
 - People below the age of 10 or above the age of 80.
 - People with athletic or muscular body types (body builders).
 - Pregnant women, people with oedema.
- For the above groups, measurements are for reference only.
- Please place the scale on a clean and levelled surface.
- Please do not use the scale in extreme high or low temperature
- Please do not jump or stamp on the scale and avoid dropping them, as the scale is a high precision electronic instrument.

- To prevent scratches on the glass panel we advise footwear not be worn while measurements are taken.
- When cleaning the scale, please use a cloth only and do not use any chemical cleaners to avoid any damage to the scale.
- This scale is designed for domestic use only and should not be used commercially for any reason at all. The (Body Fat Ratio Chart) is for reference purposes only.
- Please remove the batteries from the scale if they are not used frequently.

2. OPERATING INSTRUCTIONS

2.1. Product Description

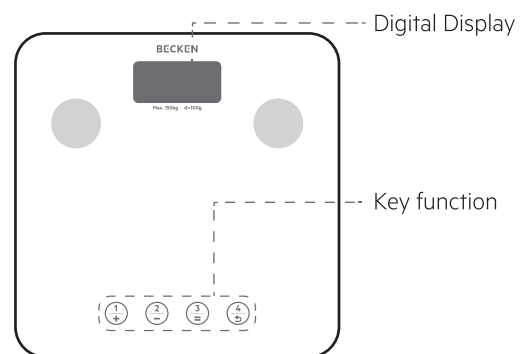Touch Keys

Press ① to increase a parameter or change the gender.

Press ② to reduce a parameter or change the gender.

Press ③ to confirm the input data.

Press ④ to go back to last step in case of any mistake in any of the operations above. Press and hold any of the buttons if you want to adjust the parameter quickly.

[you can trust]

BECKEN

2.2. Using the Product “One-key operation” – fast measurement

Select a number key from ①/②/③/④ and input your personal data into the memory before using the scale for the first time. From the second time you use the scale, you can press ①/②/③/④ for few seconds to begin measuring directly.

How to use the scale

Before using the scale, place it on a hard-flat surface. Before stepping on, there are two modes to select from.

MODE 1: Step on mode

For weighing only, no composition analysis.
Just step on the scale directly, it will automatically display your weight.


MODE 2: Switch on mode

For weighing and composition analysis.
Step 1: Switching on the scale and Family user selection.
Press and hold any of the keys ① ② ③ ④ to switch on the scale and simultaneously select the user by number.


▲ NOTE: When “0.0” appears, step on the scale with clean and dry bare feet. You must place your feet evenly on the electrodes.

Step 2: Unit setting.
When the key is flashing,
Press ① or ② to select the desired units from one of the following: kg, lb or st:lb.
Press ③ to confirm

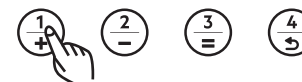Step 3: Gender selection.
When ♀ is flashing
Press ① or ② to select male or female
Press ③ to confirm (♂=Male ♀=Female)

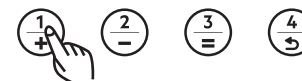Step 4: Height setting.
When the number representing the user's height is flashing
Press ① or ② to change the user's height
Press ③ to confirm
Press Unit key to change between cm/ft.

BECKEN

Step 5: Age setting.


When the number representing the user's age is blinking

Press ① or ② to change the user's age

Press ③ to confirm

Press and hold any of the buttons you want to adjust the parameter quickly.

Step 6: Weighing and measuring.

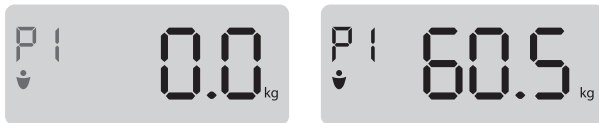When "0..." is displayed from left to right, composition analysis is underway.

Please do not step off the scale until your weight is locked on the screen.

Body fat result: The emoticons displayed below let you know your body fat result.

Underfat	Healthy
Overweight	Overfat

Reference table

Body fat percentage (%)

Age	Female				Male			
	Underfat (-)	Healthy (0)	Overfat (+)	Obese (++)	Underfat (-)	Healthy (0)	Overfat (+)	Obese (++)
10-20	≤19.5%	19.6-33.5%	33.6-38.5%	≥38.6%	≤9.5%	9.6-20.5%	20.6-25.5%	≥25.6%
21-30	≤20.0%	20.1-34.0%	34.1-39.0%	≥39.1%	≤10.0%	10.1-21.0%	21.1-26.0%	≥26.1%
31-40	≤20.5%	20.6-34.5%	34.6-39.5%	≥39.6%	≤10.5%	10.6-21.5%	21.6-26.5%	≥26.6%
41-50	≤21.0%	21.1-35.0%	35.1-40.0%	≥40.1%	≤11.0%	11.1-22.0%	22.1-27.0%	≥27.1%
51-60	≤21.5%	21.6-35.5%	35.6-40.5%	≥40.6%	≤12.0%	12.1-23.0%	23.1-28.0%	≥28.1%
61-70	≤22.0%	22.1-36.0%	36.1-41.0%	≥41.1%	≤13.0%	13.1-24.0%	24.1-29.0%	≥29.1%
71-80	≤22.5%	22.6-36.5%	36.6-41.5%	≥41.6%	≤14.0%	14.1-25.0%	25.1-30.0%	≥30.1%

Body water (hydration) percentage (%)

Age	Female				Male			
	Low (-)	Normal (0)	Overhigh (+)	High (++)	Low (-)	Normal (0)	Overhigh (+)	High (++)
10-20	≤56.0%	56.1-62.0%	62.1-68.0%	≥68.1%	≤58.0%	58.1-65.0%	65.1-72.0%	≥72.1%
21-30	≤54.0%	54.1-60.0%	60.1-66.0%	≥66.1%	≤56.0%	56.1-63.0%	63.1-70.0%	≥70.1%
31-40	≤51.0%	51.1-57.0%	57.1-63.0%	≥63.1%	≤54.0%	54.1-61.0%	61.1-68.0%	≥68.1%
41-50	≤48.0%	48.1-54.0%	54.1-60.0%	≥60.1%	≤51.0%	51.1-58.0%	58.1-65.0%	≥65.1%
51-60	≤45.0%	45.1-51.0%	51.1-57.0%	≥57.1%	≤48.0%	48.1-55.0%	55.1-62.0%	≥62.1%
61-70	≤42.0%	42.1-48.0%	48.1-54.0%	≥54.1%	≤45.0%	45.1-52.0%	52.1-59.0%	≥59.1%
71-80	≤38.0%	38.1-44.0%	44.1-50.0%	≥50.1%	≤42.0%	42.1-49.0%	49.1-56.0%	≥56.1%

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BMI (body mass index) = weight (kg) / (height)2 (m2) %

BODY	Underfat	Standard	Overfat	Obese
BMI	< 16.5	< 16.5-25	25-30	> 30

Visceral Fat (V-fat)

	Standard	Higher	Highest
Viscera-fat	≤9	10-14	≥15

2.3. Troubleshooting

Why does the scale power off after being on for just a few seconds?

Check the battery. Ensure the battery is new and inserted in the correct polarity alignment. Ensure also that your feet make good contact with the metal sheet on the scale.

If “----” is displayed on the screen, then there is something wrong with the scale.

The scale is automatically calibrated to zero (auto-zero). Please wait a few seconds until “0” is displayed. No stored parameters were found, please input your personal data.


Why does the scale only display my weight but not my body fat percentage?

When measuring body fat, movement will interfere with the measurement. Please keep your body straight and still when weighing and ensure that your feet are bare and in good contact with the metal sheet on the scale.

Why is my body fat percentage much higher than normal?

Please double-check your personal parameters given your gender, height and age. Ensure your feet are bare, clear and dry, and in good contact with the metal sheet on the scale. Also, note that any large loss of water, if you are perfectly still or not, or whether you have any calluses on your feet or not will affect the measurement of body fat.

Error messages

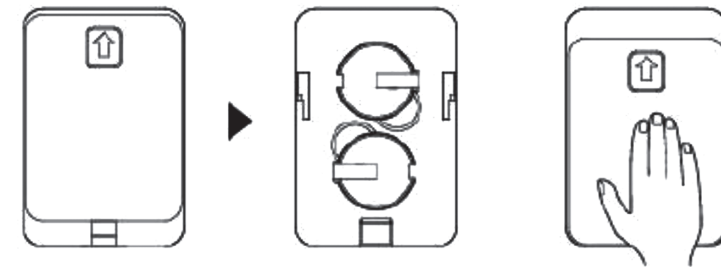Low power. Please insert the new batteries.

Weight over-load (>150 kg).

Incorrect operation caused an error in the measurement of body fat percentage.

3. TECHNICAL SPECIFICATIONS

Power supply

Battery: 2x3V CR2032 lithium battery (included)

How to replace the batteries?

When the LCD displays the “Lo” error message, turn the scale over to access the battery compartment. Slide the battery cover off, remove the old batteries and insert the new ones in the correct polarity alignment.

Specifications

The scale is equipped with a high precision strain gauge sensor system

Max. weight: 150kgs/330lbs

Min. effective weight: 8kg/17.6lb

Weight division: 0.1kg/0.2lb

Weight units: KG/LB/ST: LB

BMI/Body fat/Muscle/Water division: 0.1%

Body bone division: 0.1kg

Age range for body fat: 10-80 years

Height range: 100-220cm/3ft-3in-7ft-3in

Bone weight range: 0-15kg

Muscle range: 13.0-100.0%

Fat range: 5-50%

Water range: 35-75%

Height units: cm/ft

Memory: 4 persons

Working environment: 10-40 degree, RH<90%

Step on/Switch on

Overload/Low power indication

Auto zero resetting/Auto power off

Power: 2x3V CR2032 lithium battery (included)

4. POST-SALE SERVICE

BECKEN has designed this appliance in order to guarantee maximum reliability. However, if a problem occurs do not attempt to open the appliance - there is a risk of electrical shock. If you are not able to solve the problem, please contact the Customer Support Service of any Worten, Modelo or Continente store in order to be assisted. Worten Equipamentos para o Lar provides a 2-year warranty for manufacturing defects from the date of purchase of the appliance.

WARNING: Any attempt to repair the appliance without contacting the manufacturer or seller will void the warranty.

All users must be familiar with these instructions. Knowing what could go wrong can help prevent problems from occurring.

5. ENVIRONMENTAL PROTECTION

In order to protect the environment we strive to reduce the volume of packaging used, limiting it to three easily sorted materials: cardboard, paper and plastic. The appliance is manufactured using recyclable materials after disassembled by a specialized company. Please comply with local regulations concerning the recycling of all materials.

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