BECKEN

[you can trust]

Water Flosser

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BWF7617

Dear client,

Thank you for purchasing this product.

The quality of this appliance has been certified by several laboratory tests. This instructions manual is provided to simplify the use of the appliance as well as to guarantee a safe use.

Please read carefully the instructions and safety measures described in this manual.

We sincerely hope you will be satisfied when using your new appliance.

INDEX

1. SAFETY INSTRUCTIONS	4
1.1. Power Supply	5
1.2. Power Cord and Other Cables	5
1.3. Humidity and Water	5
1.4. Cleaning	6
1.5. General Precautions	6
2. OPERATING INSTRUCTIONS	8
2.1. Product Description	8
2.2. Product Usage	9
2.3. Maintenance	12
2.4. Troubleshooting	14
3. POST-SALE SERVICE	14
4. ENVIRONMENTAL PROTECTION	15

BECKEN

1. SAFETY INSTRUCTIONS

	WARNING	
	RISK OF ELECTRIC SHOCK DO NOT OPEN	
Warning: to avoid the risk of electric shock, do not remove the cover. The interior does not contain any parts that can be repaired by the user. Ask for assistance from the Customer Support Service.		

Important Warnings

Before using your appliance read these instructions manual carefully and keep it in a safe place as you may need to consult it in the future.

Before turning on your appliance, and to ensure a safe and correct operation, do not use if:

- It has fallen to the ground;
- You notice any evidence of damage;
- Any anomaly occurs during its operation.

To avoid hazards or deterioration of your appliance through incorrect use it is advisable to read these instructions thoroughly.

Your appliance should not be used for other purposes than the ones intended and solely for household use. Any damage resulting from using the appliance outside this scope, any incorrect use, as well as any manual modifications to the product will not be covered and automatically voids the warranty.

Assistance

In case of malfunction, take your appliance to the Customer Support Service for assistance.

This appliance can be used by children aged from 8 years and above if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised. Keep the appliance and its cord out of reach of children aged less than 8 years.

Appliances can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance.

These instructions are also available at **www.suporteworten.pt**

1.1. Power Supply

Make sure the main voltage is compatible with the appliance's. If this is not the case, please consult the Customer Support Service.

1.2. Power Cord and Other Cables

Do not tangle the cord or pull it to disconnect the appliance. Also, keep the cord away from warm surfaces.

Never use the appliance if it has a damaged cable or socket.

Take it to the Customer Support Service for replacement.

1.3. Humidity and Water

Do not place your appliance under water (e.g. for cleaning).

BECKEN

1.4. Cleaning

Before cleaning, be sure to switch off the appliance and remove the plug from the main power supply.

To clean the exterior of the product, use a dry soft cloth only.

1.5. General Precautions

- Please follow the instructions in the user manual which will help you to properly install and operate this product and enjoy its advanced features. Please keep the instructions for future reference.
- Maintenance and service operations are only to be carried out by authorized professionals.
- Any damage caused by manual modifications to the product is not covered by the warranty.
- If this product is used for other purposes than the ones originally intended, or if it is not correctly connected, the product may become damaged and the warranty will be voided.
- Never use this appliance unless you understand the operation process and the safety description to avoid being hurt due to mishandle of the appliance.
- Do not dismantle this appliance or change its components to avoid this appliance being damaged or causing electric shock.
- Do not fill the water tank with water exceeding 40 °C for oral cleaning to avoid being scalded.
- Please unplug and then remove the appliance from the outlet if it drops into water when being charged.
- Please replace the charger and then operate it again if the power cord or plug is damaged.


- Do not charge the appliance or use of it at high-temperature places.
- Do not press gums or teeth by the nozzle with excessive force to avoid being hurt.
- Put the appliance in a safe place to avoid it being damaged or dropped into water.
- Do not immerse the appliance into water to avoid malfunction or causing fire, electric shock and explosion.
- Do not plug or unplug the appliance with wet hands.
- Do not operate any appliance with a damaged charger, cord or plug, after the appliance malfunctions, if it has been damaged in any manner or if it has been dropped into water. Return the appliance to the Customer Support Service for examination, repair or adjustment.
- Do not damage, destroy or modify the power cord. Do not drag, bend or warp the power cord, and do not put it under heavy objects or clamp it among objects.
- When unplugging the appliance, grab the plug by the power cord; do not drag the cord.
- Keep the appliance and its power cord away from hot objects.
- Unplug the appliance before performing maintenance.
- Do not put the appliance in a high-temperature or moisture environment and avoid placing it exposed to direct sunlight.
- Do not charge the appliance in a high-temperature or moisture environment such as bathroom.
- Unplug the appliance if not used for a long time to avoid fire or danger.

BECKEN

- Please do not use the appliance without the nozzle attached.
- Do not load the appliance with the mouth rinse, as it may corrode the appliance.
- Please consult a medical professional to know if you can use the appliance if you've have had dentistry surgery in past two months or are receiving dental treatment.

2. OPERATING INSTRUCTIONS

2.1. Product Description

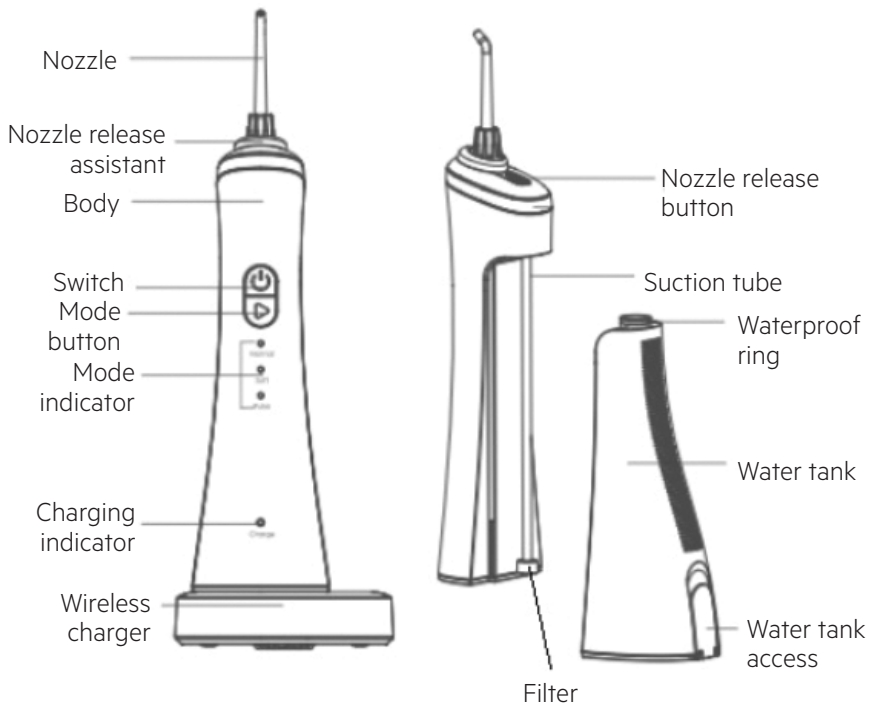
2.2. Product Usage

OPERATION INSTRUCTION

1. Check and make sure that your hands are dry before operating the appliance. Then, connect it to a power source.

2. Place the body onto the charger. The charging indicator will blink and continue blinking during charging. When the water flosser is fully charged, the charging indicator will be permanently lit. It takes about 24 hours for the water flosser to fully charge.

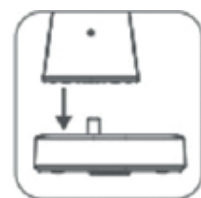NOTE: The water flosser needs to be charged if the charging indicator blinks before or after operating.

3. Remove the water flosser from the charger and then insert the nozzle into the body. Press the nozzle release button to remove the nozzle.

4. Add water to the water tank by following the two ways as below:

- 4.1. Open the water tank door and hold the water tank horizontally, add water and then close the water tank door.

BECKEN

- 4.2. Remove the water tank and hold it vertically, add water and attach the water tank to the body.

5. Select your required mode using the 'Mode button'. There are three modes to choose from, including '**Normal**', '**Soft**' and '**Pulse**' (see below). You can press the 'Mode button' to switch between the operating modes.

- A. **Normal:** For a normal teeth cleaning.
- B. **Soft:** For sensitive teeth, when the user uses the water flosser for the first time or only needs a soft cleaning.
- C. **Pulse:** For massaging gums to promote blood circulation and strengthen the teeth.

NOTE: When using the water flosser for the first time, choose the soft mode for a week. Once you have adapted to it, use in any other mode you like.

6. Hold the water flosser vertically and insert the nozzle into your mouth, aiming simultaneously at your gums and teeth, slightly opening your mouth to let the water flow out.

7. Turn on the switch to activate the water flosser. Pay attention to the

[you can trust]

following steps during operation.

- A. The water must flow out in a vertical direction towards the gums or teeth.
- B. Slowly move the water flosser along the teeth and clean them in order. Do not spray water along the periodontal pocket directly.

- 8. After operation, press the switch to turn off the appliance. The appliance will shut off automatically after operating for 2 minutes. Please press the switch again to continue operating it.
- 9. Once you have finished, you can follow the operation as below to remove the residual water:

- A. Open the water tank door and empty the water tank.

- B. Press the switch; the water flosser will start working to automatically drain out the residual water inside.

- C. Press the switch again to turn off the water flosser and then close the water tank door.

BECKEN

D. Wipe the water flosser with a dry cloth.

NOTE: The water flosser must be free of residual water to avoid bacterial infection. Therefore, if it is not used for a long time, please clean the water flosser well, then dry and store it.

2.3. Maintenance

Clean the water flosser with clean water or neutral detergent, but never use corrosive or abrasive detergent as it will damage the water flosser. Do not wash it with water hotter than 50 °C.

1. Clean the body

- A. Clean the body with a dry or wet cloth.
- B. Do not immerse the body into water for cleaning.

2. Clean the water tank

- A. Wipe the water tank with a dry or wet cloth, or detach it and rinse it with water.
- B. Empty the water tank completely if it is not used for a long time.

3. Clean the nozzle

- A. Wipe the nozzle with a dry or wet cloth, or detach it and rinse it with water.
- B. Do not bend, pull and warp the nozzle.

4. Clean the charger

Wipe the metal part of the plug with a dry cloth every 6 months to avoid dust from accumulating and resulting in a poor contact.

DETACH RECHARGEABLE BATTERY

Do not dismantle the water flosser unless it needs to be discarded, or it

will be damaged. Unplug the appliance before removing the battery.

1. Make sure the battery has been exhausted before removing it. To run out the battery, press the switch to activate the water flosser and let it run until it stops working automatically. Repeat the operation for several times until the battery dies.

2. Remove the water tank and unscrew the two screws at the bottom cover with a screwdriver.

3. Remove the two screws at the top cover with a screwdriver, and then detach the top cover and remove the unit core.

4. Unscrew the screws at the bracket of the unit's core with a screwdriver.

5. Shear off the connecting wire of battery with scissors, and then remove the battery.

NOTE: Do not short-circuit the battery and never energise the water flosser again after removing its battery.

BECKEN

2.4. Troubleshooting

Problem	Possible cause	Solution
The water flosser does not work.	The water flosser is being used for the first time or it has not been used for more than 3 months.	Charge the water flosser for 24 hours.
Water pressure is low when working.	There is something blocking the nozzle. The nozzle is deformed	Replace the nozzle.
	The filter is blocked.	Clean the filter.
The water flosser can only be operated for a short time after being fully charged.	The battery has reached its end life.	Contact the Customer Support Service.
	The water flosser has not been properly charged.	Charge the water flosser for more than 24 hours before operation.
The charging indicator does not illuminate when charging.	The power plug has not been properly connected.	Connect the power source properly.
	Circuit malfunction.	Contact the Customer Support Service.

3. POST-SALE SERVICE

Becken has designed this appliance in order to guarantee maximum reliability. However, if a problem occurs do not attempt to open the appliance - there is a risk of electrical shock. If you are not able to solve the problem, please contact the Customer Support Service of any Worten, Modelo or Continente store in order to be assisted. Worten Equipamentos para o Lar provides a warranty against manufacturing defects for a 3-year period, or for another period provided for by the applicable legislation, from the date of purchase of the appliance.

WARNING: any attempt to repair the appliance without contacting the manufacturer or seller will void the warranty.

All users must be familiar with these instructions. Knowing what could go wrong can help prevent problems from occurring.

4. ENVIRONMENTAL PROTECTION

In order to protect the environment, we strive to reduce the volume of packaging used, limiting it to three easily sorted materials: cardboard, paper and plastic. The appliance is manufactured using recyclable materials after disassembled by a specialized company. Please comply with local regulations concerning the recycling of all materials.

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Worten Equipamentos p/ o Lar S.A.